# DIENTES DE RESINA ACRÍLICA

## USO PREVISTO

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# ACRYLIC RESIN TEETH

## EXPECTED USE

Acrylic resin teeth are made of polymethyl methacrylate and are medical devices designed to replace the loss of natural teeth. They are used in the dental sector for the manufacture of dental prostheses, the most common being the muco-supported (total removable prosthesis), and the dento-muco-supported prosthesis (removable partial denture of acrylic resin and removable partial metal prosthesis). In some cases they are used for implant-supported prostheses (over dentures or hybrid prostheses). They are intended for all those who have suffered loss of natural teeth due to injuries such as cavities, gum disease or any other health problem, in order to recover the functionality and well-being of the patient.

## COMPOSITION

Polymethyl methacrylate  
Ethylene glycol dimethacrylate (cross-link)  
Fluorescence  
Pigments

## FEATURES

Acrylic resin teeth are classified into three categories according to their structure: two, three and four layers, likewise two classes are identified: class I that corresponds to the anterior teeth that go from canine to canine, and class II that go from the first premolar to the second molar in both hemiarcs.

The product has a wide variety of colors and shapes in both anterior and posterior teeth that can be adapted to the clinical condition and physiognomy of each patient.

Acrylic resin anterior teeth are grouped into three basic shapes: square, ovoid, and triangular according to the patient's facial contour. 2 and 3-layer posterior teeth can vary in their occlusal surfaces, inverted cusp molars (0 degrees), semi-anatomical molars (10 and 20 degrees), anatomical molars (33 degrees), and 4-layer molars can vary by type of occlusion and morphological characteristics that adapt to a normal and crossed joint.

## INSTRUCTIONS FOR USE

The dentist is the one who decides the color, size and type of joint according to the physiognomy and clinical conditions of the patient. This is done based on the mold chart and shade guide provided by New Stetic.

### Selection of anterior teeth

Select teeth directly on the patient to suit their anatomical and facial characteristics.

- Square facial contour: Teeth where the angular lines of the crown are practically parallel, the incisal edge is straight and is accompanied by a gingiva with little distal inflection. The mesio incisal angle is close to 90°, while the distal angle is somewhat more obtuse.

- Triangular facial contour: Teeth where the angular lines converge sharply towards the gingival. The incisal edge is straight and the mesioincisal and disto incisal angles are rounded. The buccal surface is slightly convex and the incisal third is flattened.

- Oval facial contour: Teeth where both the angular lines and the buccal surface are very convex.

### Selection of posterior teeth

A variety of occlusal sizes and shapes are stated that give rise to cusp tilt, this is one of many factors in the laws of articulation to achieve balanced occlusion.

To select the cusp angulation of the posterior teeth in the different edentulous patients, the biological, anatomical and physical characteristics must be analyzed and related, such as: Alveolar ridge resorption, TMJ limitations, intermaxillary relationship (class I, II, III and cross bites), neuromuscular coordination, among others.

- Anatomical posterior teeth with 33° cusp angulations are the most favorable for balanced occlusion, artificial teeth with this inclination are comparable with natural teeth, they have very marked cusps and fossae that interrelate correctly.

- The 20° teeth are semi-anatomical, they have a cusp angulation that provides less height to make balanced contacts in eccentric excursions of the jaw than the 33° anatomical tooth.

- 0° posterior teeth or non-anatomical teeth have a flat occlusal face, without cusps, but with some depressions that allow a certain chewing efficiency. Its application is effective when it is difficult or impossible to accurately record the centric relationship of the patient or when there are abnormal mandibular relationships in patients with little ability to adapt to dental prostheses.

**Note:** Often the surfaces of artificial teeth must be modified to fit the balanced occlusal scheme of each patient.

Color selection should be done on the patient under appropriate lighting conditions.

The selection of colors is made with the help of the New Stetic A-D shade guide according to each brand.

The dental laboratory designs the prosthetic structure based on the dentist's indications and the study model, following the instructions for use of the other materials required to obtain the final device.

The transfer of patient-specific intermaxillary relationships must be translated into a semi-adjustable articulator, where individual patient values can be included, it is an essential step to perform functionally effective dental prostheses.

The dentist and dental laboratory must ensure a balanced occlusion, to achieve this it is important that basic concepts of occlusion are handled, not taking them into account implies altering the harmonic balance that must exist between the components of the stomatognathic system, the integrity and useful life of the dental prostheses.

## WORKING TOGETHER WITH OTHER DEVICES

This product is used in combination with acrylic resins for dental prostheses. In some types of restorations, the teeth work with noble metal frameworks.

## RESIDUAL RISKS ASSOCIATED WITH THE USE OF THE PRODUCT

The chemical properties of the material ensure adequate chemical adhesion and allow compatibility between the tooth and the acrylic resin for the prosthetic base. In some types of restorations, such as skeletal partial dentures, the teeth work together with metal structures, but the retention between the two devices is mechanical, avoiding in both cases the separation of the tooth from the base of the prosthesis.

Artificial teeth have specific characteristics, shapes and sizes to adapt to the clinical needs of each patient, making inappropriate wear unnecessary and facilitating the planning of a balanced occlusion by the dentist and dental laboratory. Excessive wear and inadequate clinical planning can affect the functionality of the prosthesis causing fracture of the teeth, so the prosthesis must be replaced or repaired.

The physical properties of the material allow a smooth surface to be recovered using conventional polishing techniques, thus avoiding plaque build-up. In addition, the patient must follow the recommendations given by the dentist to establish daily cleaning routines.

## CONTRAINDICATIONS

Do not use in case of proven allergy to polymethyl methacrylate.

## WARNINGS

Wear a face mask or mask and, if possible, suction equipment to avoid breathing dust generated during drilling.

## PRECAUTIONS

Handling of the device must be in accordance with the principles of good hygiene practice.

Avoid making excessive adjustments during milling to ensure the functionality and shade effect of the teeth. Treated tooth surfaces must show a polished and shiny surface finish prior to mouth fitting. Use silicone polishers for the final polish, to avoid accumulation of bacterial plaque.

The acrylic teeth must not come into contact with direct heat or a flame (eg from an alcohol lighter), this removes the shine from the surface and prevents it from recovering it.

The dental prosthesis must be cleaned with a conventional product. Avoid contact with solvents.

The use of toothpaste to clean the prosthesis is contraindicated, since it causes the opacity and deterioration of the outer layer of the devices in the long term, promoting the accumulation of bacterial plaque and tartar.

## RECOMMENDATIONS

For best results, the use of acrylic resins manufactured by New Stetic is recommended.

## STORAGE

The product should be stored in a cool and dry place, at a temperature no higher than 30 °C.

## FINAL DISPOSITION

Do not dump into bodies of water. The final disposal of the product must be done in accordance with current local regulations.

## USEFULL LIFETIME.

10 years

Health Certificate: Colombia

INVIMA2015DM-001484-R2

## NORMATIVITY

Acrylic Resin Teeth product complies with ISO 22112 standard


In case of serious accidents, contact: dosorio@newstetic.com (Technical Director of Medical Devices) and the email infocolombia@newstetic.com.

## NORMATIVE DIMENSIONS

A Length of central without collar.


B Width of central.

C Width of 6 anteriors on a straight line

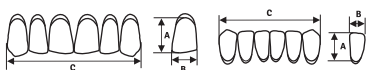A-B Width of 4 upper or lower posteriors on a straight line.

C-D Width of the upper first molar or lower left molar
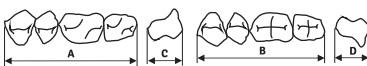
## SYMBOL GLOSSARY

|   |                      |
|---|----------------------|
| Keep away from sunlight or direct heat sources. | Keep dry             |
| Maximum storage temperature.                    | Instructions for use |

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**New Stetic S.A.**